HS

ES|EN|FR

BALANZA DE PRECISIÓN
PRECISION SCALE
BALANCE DE PRÉCISION

V.2.3
09/02/2024

ES

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EN

1. TECHNICAL DATA	28
2. DESCRIPTION OF THE DEVICE	29
2.1 COMPONENTS	29
2.2 DESCRIPTION OF THE KEYPAD	29
2.3 POSSIBLE INDICATIONS	30
3. BASIC INFORMATION (GENERAL INFORMATION)	30
3.1 INTENDED USE	30
3.2 IMPROPER USE	30
3.3 GUARANTEE	31
3.4 MONITORING THE MEANS OF CONTROL	31
4. BASIC SAFETY RECOMMENDATIONS	31
4.1 FOLLOW THE RECOMMENDATIONS IN THE OPERATING INSTRUCTIONS	31
4.2 STAFF TRAINING	31
5. TRANSPORT AND STORAGE	31
5.1 CHECKING UPON RECEPTION	31
5.2 PACKAGING/RETURN	31
6. UNPACKING, LOCATION AND START-UP	32
6.1 LOCATION AND PLACE OF OPERATION	32
6.2 UNPACKING, CHECKING AND LOCATION	32
6.3 LEVELLING	33
6.4 ELECTRICAL SUPPLY	33
6.5 FIRST START-UP	34
6.6 CONNECTION OF AUXILIARY DEVICES	34
7. SETTINGS	34
7.1 MODELS HS120I - HS210I	34
7.1.1 AUTOMATIC ADJUSTMENT USING INTERNAL WEIGHT	34
7.1.2 AUTOMATIC ADJUSTMENT IS CARRIED OUT	34
7.1.3 MANUAL START OF THE INTERNAL SETTINGS	35
7.2 MODELS HS120E - HS210E	36

8. BASIC MODE	37
8.1 SWITCHING ON THE DEVICE	37
8.2 SWITCHING TO STAND-BY MODE	37
8.3 ZERO SETTING	38
8.4 SIMPLE WEIGHING	38
8.5 CHANGING THE UNIT	38
8.6 TARE WEIGHING	38
9. APPLICATIONS	39
9.1 PARTS COUNT	39
9.2 DETERMINING THE PERCENTAGE	40
10.MENU	41
10.1 MENU NAVIGATION	41
10.2 DESCRIPTION OF THE MENU	43
11.RS-232C INTERFACE	44
11.1 TECHNICAL DATA	44
11.2 INTERFACE CABLE	44
11.3 INTERFACE PARAMETERS	44
11.3.1 BAUD RATE	44
11.3.2 CONDITION OF EDITING/PRINTING DATA	45
11.4 REMOTE CONTROL COMMANDS	45
11.5 FRAME IN STABLE MODE, STABLE EVERY 3S AND MANUAL	45
12.MAINTENANCE	46
12.1 CLEANING	46
12.2 MAINTENANCE, CONSERVATION IN GOOD WORKING ORDER	46
12.3 WASTE TREATMENT	46
13 HELP IN CASE OF MINOR BREAKDOWNS	46
13.1 ERROR MESSAGES	47

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1. TECHNICAL DATA

MODEL	HS120E		HS210E	HS120I	HS210I
Default unit	g				
Minimum graduation (d)	0.0001 g				
Weighing range (Max.)	120 g	210 g		120 g	210 g
Reproducibility	0.0002 g				
Linearity	± 0.0004 g				
Calibration weight	100 g (E2)	200 g (E2)		internal	
	recommended (class) not included in the delivery				
Possible adjustment points	100 g	100 g		-	
	-	200 g		-	
Growth time of the signal	3 s			4 s	
Preparation time	8 h				
Weighing units	g, oz, ct, lb				
Minimum unit mass per component for the unit count	1 mg (under laboratory conditions)				
	10 mg (under normal conditions)				
Number of reference components during the unit count	10, 20, 50,1000			10, 20, 50,100, 1000	
Stainless steel weighing plate	Ø 90 mm				
Housing dimensions (WxDxH) [mm]	230 x 310 x 330				
Net weight [kg]	4.4			4.9	
Environmental conditions	from +10°C to +30°C				
Relative humidity in the air	20~85% (without condensation)				
Mains adapter input voltage	AC100 - 240V, 50 - 60Hz 2.0A				
Input voltage scale	DC 12 V, 2 A				
Interface	RS-232C				

To avoid unnecessary damage to the scale, please pay close attention to the following recommendations:

- Do not use the scale in hazardous areas.
- Disconnect the scale from the power supply before connecting it to (or disconnecting it from) any other installation.
- Excessive electromagnetic interference may cause the display to read incorrectly. Once any electromagnetic interference is reduced, the scale should operate normally.
- Avoid spilling any liquid on the surface of the scale. Use a slightly damp cloth to clean the scale

2. DESCRIPTION OF THE DEVICE

2.1 COMPONENTS

NO.	NAME
1	GLASS WINDSHIELD
2	RS232 INTERFACE
3	WEIGHING PLATE
4	BUBBLE LEVEL
5	SCREEN
6	KEYPAD
7	HEIGHT-ADJUSTABLE LEG
8	WINDSHIELD HANDLE
9	POWER SUPPLY PLUG

2.2 DESCRIPTION OF THE KEYPAD

ON/OFF key
On/Off
Exit the menu

MODE key
Switching between weighing units.
Select the application

CAL key
Adjustment (calibration)

MENU key
Saving the reference value


TARE key
Tare
Set to zero
Saving settings

PRINT key
Changing the menu settings
Printing the edited value

2.3 POSSIBLE INDICATIONS

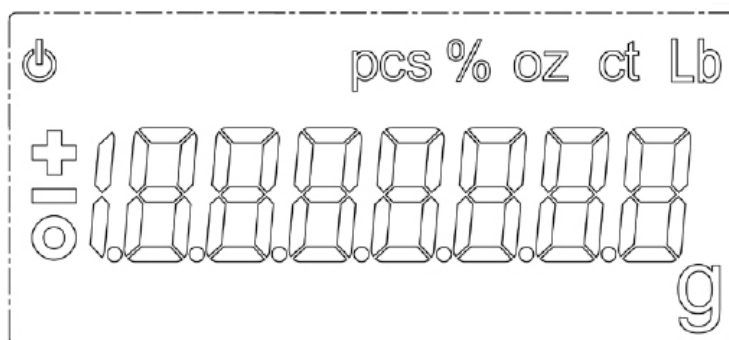INDICATION	DESCRIPTION
0	stable value indicator
	scale in standby mode (STAND-BY)
g	Unit of weight "gram"
ct	Unit of weight "carat"
lb	Unit of weight "pound"
oz	Unit of weight "ounce"
pcs	Unit counting application
%	Percentage determination application

3. BASIC INFORMATION (GENERAL INFORMATION)

3.1 INTENDED USE

The scale you have just purchased is used to define the mass (the weighing value) of the weighed material. It should be treated as a "non-automatic" scale, i.e. the material to be weighed must be carefully placed manually in the centre of the plate. The mass value is read after the scale has been stabilized.

3.2 IMPROPER USE

This scale is not intended for dynamic weighing, meaning even small amounts of the heavy material are added or removed during use. Due to the existence of the "stabilization compensation" mechanism, the scale may indicate incorrect values! (Example: slow loss of liquid from the container placed on the scale).

Do not load the weighing plate for an extended period of time. Otherwise, the measuring mechanism may become damaged.

Avoid any knocks or overloading of the plate above the maximum load (Max.), including the load used in the tare. Otherwise, the scale may be damaged.

Never use the scale in a place where there is a risk of explosion. The standard version has no explosion protection.

No structural changes should be made to the scale. Any attempt to change the scale may result in errors in the display and control unit, violations of safety regulations or damage to the equipment.

The balance may only be used in accordance with the recommendations described.

3.3 GUARANTEE

The guarantee is void in case of:

- not following the recommendations in the instruction manual;
- not using the device in accordance with the applications described;
- modifying or opening up the device;
- mechanically damage, or damage the device through the action of supplies, liquids; normal wear and tear;
- placing the device improperly or using the device with inappropriate electrical installation;
- overloading the measuring mechanism.

3.4 MONITORING THE MEANS OF CONTROL

Within the framework of the quality control system, the technical measurement properties of the scale and, if accessible, of the test weight must be checked regularly. To this end, the user in charge must define the appropriate frequency as well as the standard and limits of these controls.

4. BASIC SAFETY RECOMMENDATIONS

4.1 FOLLOW THE RECOMMENDATIONS IN THE OPERATING INSTRUCTIONS

Please read these operating instructions carefully before installing and operating the scale, even if you have used BAXTRAN scales before.

4.2 STAFF TRAINING

This device may only be used and maintained by trained personnel.

5. TRANSPORT AND STORAGE

5.1 CHECKING UPON RECEPTION

Immediately after receipt of the shipment, it is essential to check that the packaging is not visibly damaged. The same procedure applies to the device after it has been removed from its packaging.

5.2 PACKAGING/RETURN

This device may only be used and maintained by trained personnel.

- All components of the original packaging must be kept in case the device has to be returned.
- The transport for returning the device must always be in the original packaging.
- Before sending the device, all connected cables and the individual/mobile units must be disconnected.
- Unpack the plate and accessories separately, as well as the power adapter.
- Secure the glass protection housing against slipping (e.g. by using adhesive tape).
- All parts must be secured against movement and damage, as shown in the picture.

6. UNPACKING, LOCATION AND START-UP

6.1 LOCATION AND PLACE OF OPERATION

The scales are built to provide reliable measurement results under normal operating conditions.

Choose a suitable location for the scale to ensure that it works accurately and fast.

THE PRINCIPLES SET OUT BELOW MUST BE FOLLOWED AT THE PLACE OF OPERATION:

- Position the scale on a stable, flat surface.
- Avoid extremes of temperature as well as temperature changes due to, e.g., the presence of radiators or working in an area with a risk of direct exposure to sunlight.
- Protect the scale from direct drafts caused by open doors and windows.
- Avoid shaking during weighing.
- Protect the scale from high atmospheric humidity, vapour and dust.
- Do not expose the device to high humidity for a long period of time. The device may become covered with dew (condensation of atmospheric humidity) if it is moved from a cold environment to a warmer one. If this happens, the device must remain off for approximately 2 hours at room temperature in order for it to acclimatize.
- Avoid static charges between the material to be weighed and the scale's container.

In the event of electromagnetic fields, static loads or an unstable power supply, the weight display may deviate (resulting in incorrect weighing). Therefore, change the location of the scale.

6.2 UNPACKING, CHECKING AND LOCATION

Open the packaging and carefully remove all the components.

Check that all the delivery components are there and check their integrity.


DELIVERED COMPONENTS / STANDARD ACCESSORIES

- Scale
- Mains adapter
- Instruction manual

KEEP THE DEVICE DISCONNECTED FROM THE MAINS DURING INSTALLATION

Place the scale in the intended place of use. Correct positioning is essential for accurate weighing with high-resolution precision scales (see Chapter 7.1).

Remove the transport protections (models HS120I and HS210I only)

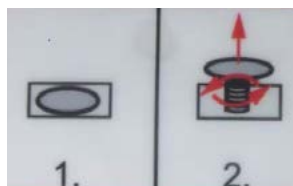
Install the following components:

- weighing plate support,
- weighing plate.

6.3 LEVELLING

Place the scale in a horizontal position using the adjustable feet with screws. The air bubble in the level must be centred.

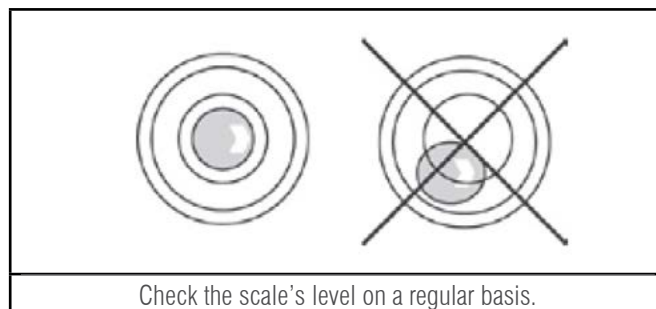6.4 ELECTRICAL SUPPLY

Select the appropriate plug for the country of use.

Check that the supply voltage is correctly selected. The device can only be connected to the power supply if the data indicated on the device (on the sticker) are identical to the local voltage.


Only use original BAXTRAN network adapters. Use of another product requires authorization from BAXTRAN.

Important:

- Check the integrity of the electric cable before using it.
- The mains adapter must not come into contact with liquids.
- The electrical outlet must be easily accessible at all times.

Connect the mains adapter to the connection socket at the rear of the scale and connect it to the power supply.

The panel will light up. As soon as the power adapter is connected and the scale is plugged in, the power symbol will appear .

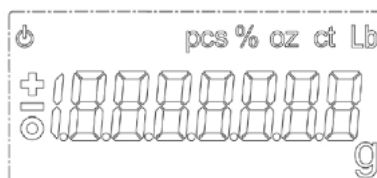i

The error message <P FAIL> is displayed when the scale has been disconnected from the power supply without pressing the button .

Solution:

Press the  key. The control panel will run the test. The scale is ready for weighing once the mass display appears.

6.5 FIRST START-UP

To ensure that the electronic scales display the correct results, they must be set to the correct operating temperature (see "Preparation Time" in Chapter 1). During preparation, the scale must be connected to the power supply (mains or battery).

The accuracy of the device depends on the gravitational acceleration.

The instructions in the "Settings" chapter must be followed.

6.6 CONNECTION OF AUXILIARY DEVICES

Before plugging or unplugging auxiliary devices (printer, computer, etc.) to/from the interface, the scale must be unplugged from the mains.

The scale may only be operated with accessories and auxiliary devices from BAXTRAN that are correctly synchronized with the scale.

7. SETTINGS

Since the value of the gravitational acceleration is not the same at all points on the Earth, each scale must be adjusted – according to the weighing principle resulting from physical principles – to the gravitational acceleration at the location where the scale is located (only if the scale has not been adjusted to the scale's location at the factory). This settings' process must be performed during the first start-up and every time the scale's location changes, as well as when the room temperature changes. To ensure accurate weighing results, we also recommend that the scale be systematically adjusted in weighing mode.

Ensure that environmental conditions are stable. Follow the required preparation time (see Chapter 1) to stabilize the scale. Make sure that the scale plate is free of objects.

7.1 MODELS HS120I - HS210I


7.1.1 AUTOMATIC ADJUSTMENT USING INTERNAL WEIGHT

Mandatory setting in the <C1-0> menu, see chapter 11.2

7.1.2 AUTOMATIC ADJUSTMENT IS CARRIED OUT

- If there has been a 2°C change in temperature since the last setting,
- After the defined time has elapsed, the option to choose between 2, 3, 4, 6, 8 hours (C7-0-C7-4, see Chapter 11.2).

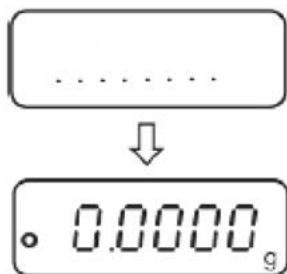In the event of a settings' error (e.g., the weighing plate is not free of objects), the display will show the error message <CAL NO>. Repeat the settings' process.

7.1.3 MANUAL START OF THE INTERNAL SETTINGS

Mandatory settings in the <C1-2> menu, see chapter 11.2

PROCEDURE

In the weighing mode press the  key

The settings are made automatically

After establishing the correct settings, the scale will automatically return to weighing mode. In the event of a settings' error (e.g., the weighing plate is not free of objects), the display will show the error message <CAL NO>. Repeat the settings' process.


7.2 MODELS HS120E - HS210E**i**

As far as possible, adjustments should be performed with a mass close to the maximum load of the scale (recommended settings' weight, see Chapter 3.1.1). 1). Information on adjustment weights is available on the website www.baxtran.com

- Ensure that environmental conditions are stable. In order to stabilize the scale, the necessary preparation time must be provided.
- Make sure that the scale plate is free of objects.
- The settings in the <C1-1> menu should be borne in mind, see chapter 11.2.

PROCEDURE


In the weighing mode press the  key.

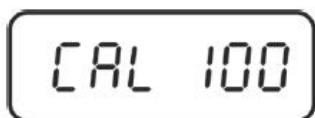The required mass of the calibration weight is displayed (e.g., 100 g).


Press the  key again, the indication <CAL in> will appear.

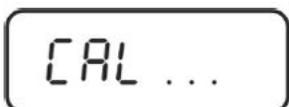
Make sure that the scale plate is free of objects.

Wait for the <CAL dn> display and then carefully place the calibration weight in the centre of the plate. Close the door of the windshield.

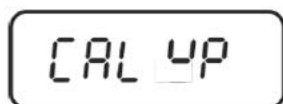After successful calibration, the indication <CAL up> will appear.

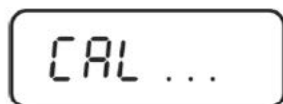Remove the calibration weight. The scale returns automatically to the weighing mode.

8. BASIC MODE

8.1 SWITCHING ON THE DEVICE

As soon as the power adapter is connected and the scale is plugged in, the power symbol will appear ⏻.

To turn on the device, press the  key

The device performs the control panel diagnosis.
The scale is ready for weighing once the mass display appears.

8.2 SWITCHING TO STAND-BY MODE

Press the  key, the display will switch off.

The power symbol appears ⏻.

- In standby mode, the scale is ready for use after it has been switched on, without any preparation time.
- To switch off the device permanently, disconnect it from the power supply.
- The scale will start up in the same mode in which it was turned off.

8.3 ZERO SETTING

To obtain optimum weighing results, the scale must be zeroed before weighing

Unload the scale.

Press the  key. Wait for the zero indication.

8.4 SIMPLE WEIGHING

Place the material to be weighed on the plate and close the door of the protective housing.

Wait for the stabilisation indication .
Read the weighing result.

Warning of overloading

It is imperative to avoid any knocks or overloading of the device above the maximum load (Max.), including the load used in the tare. Otherwise, the scale may be damaged. The letter "E" indicates overloading above the maximum values. Unload the weighing system or decrease the initial load.

8.5 CHANGING THE UNIT

Press the MODE key several times to switch the mass display between the available weighing units and to change the application.

g ⇌ oz ⇌ ct ⇌ lb ⇌ Pcs ⇌ %


8.6 TARE WEIGHING

The mass of any container used for weighing can be tared by pressing the corresponding key, and thus the net mass of the weighed material will be displayed for subsequent weighing.

Zeroing the scale

Place the empty container on the scale. Its mass will appear.

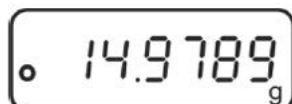Wait for the stabilization index to appear and press the key .

The display will show zero.

Weigh the material to be weighed. The indicated mass corresponds to its net mass.






- If the scale is not loaded, the stored tare value is displayed with the negative value symbol.

- To delete the stored tare value, unload the plate and press .

9. APPLICATIONS

9.1 PARTS COUNT

Before counting the units using the scale, it is necessary to define the average mass of a unit (unit mass) known as the reference value. For this it is necessary to place a certain number of units to be weighed. The scale determines the total mass and divides it by the number of units, known as number of reference units. Then, based on the average mass calculated for one part, the counts will be performed.

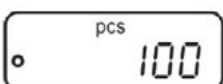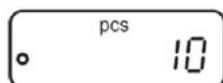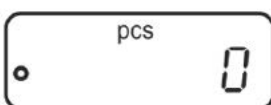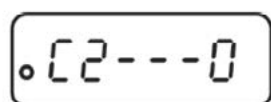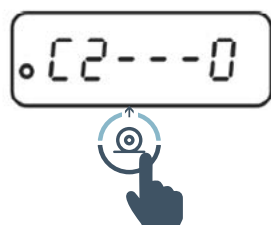CHOOSE THE NUMBER OF REFERENCE UNITS

In the weighing mode, edit menu item "C2", see Section 11.1.

Use the  key to select the desired number of reference units, see Chapter 11.2.

Validate the setting by pressing .

Exit the menu by pressing : When prompted to "SAVE", press the  key.

The scale returns automatically to the weighing mode.

EDITING THE APPLICATION

Press the  key several times until the "Pcs" symbol appears

SETTING TO ZERO / TARE

Press the  key to zero the scale, or to tare the scale container being used.

REFERENCE PARTS WEIGHING

Place the requested number of units counted according to the defined number of reference units.

Save the reference value by pressing . The scale automatically shows the average value of the part's mass. Remove the reference load. The scale is in unit-counting mode and counts the units on the plate.

PARTS COUNTING

Place the material to be weighed on the plate and read the number of parts.

SWITCHING BETWEEN UNIT AND MASS DISPLAY MODES

It is possible to switch between the display units by pressing the  key, see Chapter 8.5.

PRINTING


After connecting the optional printer, by using the key  (factory setting), the display value can be printed out.

- Take into account the minimum mass of a unit (see Chapter 1 "Technical data").
- The reference mass is also stored after the scale is switched off, until the reference value is defined again.

9.2 DETERMINING THE PERCENTAGE

Defining the percentage value allows you to see the mass in percentage form with reference to the reference mass which equals 100%.

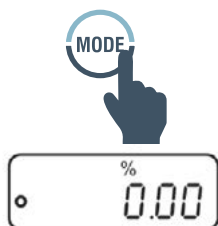EDITING THE APPLICATION

Press the  key until the "%" symbol appears.

SETTING TO ZERO / TARE


Press the  key to zero the scale, or to tare the scale container being used.

REFERENCE VALUE SETTING (100% VALUES)

Place a reference mass that corresponds to 100%.

Save the reference value by pressing . Remove the reference load.

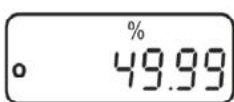DETERMINING THE PERCENTAGE

Place the material to be weighed. The mass of the sample is shown as a percentage with reference to the reference mass.

SWITCHING BETWEEN PERCENTAGE AND MASS DISPLAY MODE

It is possible to switch between the display units by pressing the  key, see Chapter 8.5.

PRINTING

After connecting the optional printer, by using the  key (factory setting), the display value can be printed out.

The reference mass (100%) is also stored after the scale is switched off, until the reference value is defined again.

10. MENU


10.1 MENU NAVIGATION

1

ENTER THE MENU

In the weighing mode press the  key, and then, .

The first menu item "C1" appears with the current setting.

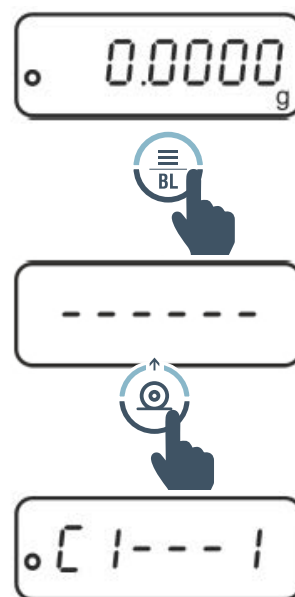2

SELECTING MENU ITEMS

Menu items with current settings can be selected in order by pressing the  key.


3

CHANGING SETTINGS

The settings in the selected menu item can be changed using the  key. After pressing the PRINT key, the following setting is displayed.



Validate using the  key. The menu item will appear. Enter the following settings or exit the menu, and save the settings (see step 4 or 5).

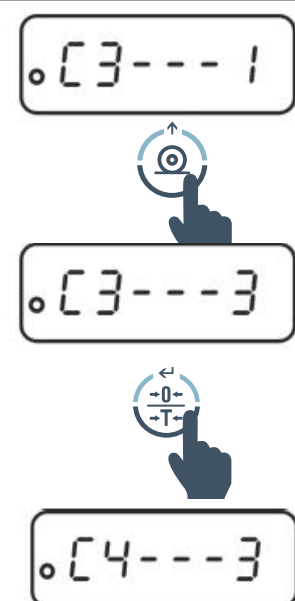
4**SAVE THE SETTINGS AND EXIT THE MENU**

Press the  key. A "SAVE" prompt will appear.

Pressing the TARE key will save the changes. The scale returns automatically to the weighing mode.

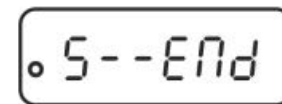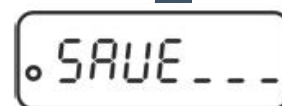**5****CANCEL**

Press the ON/OFF key again, the changes made will not be recorded. The scale returns automatically to the weighing mode.

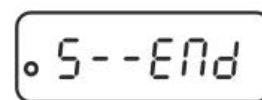10.2 DESCRIPTION OF THE MENU

MENU ITEM	SETTINGS	DESCRIPTION
Settings' mode	C1 - 0	Automatic internal setting (models HS120I and HS210I)
	C1 - 1	External setting (models HS120E and HS210E: always use this setting)
	C1 - 2	Automatic internal setting (models HS120I and HS210I)
Number of reference units	C2 - 0	10
	C2 - 1	20
	C2 - 2	50
	C2 - 3	100
	C2 - 4	1000
Automatic zero correction.	C3 - 0	The Auto Zero function is off
	C3 - 1	The range of the AutoZero function ± 1 number
	C3 - 2	The range of the AutoZero function ± 2 number
	C3 - 3	The range of the AutoZero function ± 3 number
	C3 - 4	The range of the AutoZero function ± 4 number
	C3 - 5	The range of the AutoZero function ± 5 number
	C3 - 6	Not documented
Baud rate	C4 - 0	1200
	C4 - 1	2400
	C4 - 2	4800
	C4 - 3	9600
Editing/printing data	C5 - 0	Automatic data formatting/printing of stable weight values
	C5 - 1	Using the remote control commands (P)
	C5 - 2	Continuous editing/printing of stable weight values (3-second cycle)
	C5 - 3	Editing/printing the stable and unstable weight value after pressing the PRINT key
	C5 - 4	Not documented
Keypad sound	C6 - 0	switched on
	C6 - 1*	switched off
Remaining time for automatic internal settings (models HS120I and HS210I only)	C7 - 0*	2 hours
	C7 - 1	3 hours
	C7 - 2	4 hours
	C7 - 3	6 hours
	C7 - 4	8 hours
Filter	C8 - 0*	Slow
	C8 - 1	Medium
	C8 - 2	Fast
	C8 - 3	Not documented

11. RS-232C INTERFACE


In order to connect auxiliary devices (printer, computer), the scale has the RS-232C interface as standard.

To ensure communication between the scale and auxiliary devices, the following conditions must be met:

- The scale can be connected via the interface to an auxiliary device using the appropriate cable. Only BAXTRAN's interface cables ensure error-free operation.
- The communication parameters (baud rate, bits, parity) of the scale and the auxiliary device must match.

The data is transmitted asynchronously in ASCII code.

11.1 TECHNICAL DATA

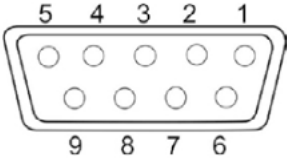PLUG	9-PIN MINIATURE CONNECTOR - D-SUB	
		Pin 2: RX (data reception) Pin 3: TX (data transmission) Pin 5: GND (mass)
BAUD RATE	possibility of choosing between 1200/2400/4800/9600	
PARITY	8 bits, no Parity, 1 stop bit, 1 start bit	270407 RS0004



11.2 INTERFACE CABLE

9-PIN SCALE CONNECTOR	9-PIN COMPUTER CONNECTOR
2	3
3	2
5	5

11.3 INTERFACE PARAMETERS

11.3.1 BAUD RATE

In this menu item, the baud rate of the data is adjusted to different auxiliary devices equipped with RS-232C. The baud rate defines the speed of data transfer over the standard interface. To ensure correct data transmission, set the scale and the auxiliary device to the same value.

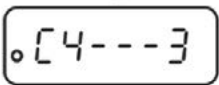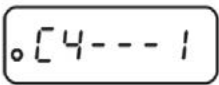Edit menu item "C4", see Chapter 10.1.

Select the desired setting by pressing


Possibility of choosing between:

- C4 - 0** 1200
- C4 - 1** 2400
- C4 - 2** 4800
- C4 - 3** 9600

Save/return to weighing mode, see Chapter 11.1

11.3.2 CONDITION OF EDITING/PRINTING DATA

In this menu item the data transmission mode is determined.

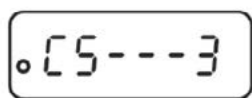Edit menu item "C5", see Chapter 10.1.

Select the desired setting by pressing

Possibility of choosing:


C5 - 0 Automatic data formatting/printing of stable weight values

C5 - 1 Using the remote control commands (P)

C5 - 2 Continuous editing/printing of stable weight values (3-second cycle)

C5 - 3 Editing/printing the stable and unstable weight value after pressing

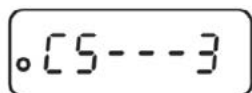C5 - 4 Not documented

Save/return to weighing mode, see Chapter 11.1

11.4 REMOTE CONTROL COMMANDS

COMMANDS		FINISHER		Function
ASCII	Hex	<CR>	<LF>	
O	4F	0D	0A	Standby ON/OFF
A	54	0D	0A	TARE
C	43	0D	0A	CAL ext
M	4D	0D	0A	Change unit
W	57	0D	0A	Printing
w	77	0D	0A	Printing

11.5 FRAME IN STABLE MODE, STABLE EVERY 3S AND MANUAL

In modes g, oz, ct or lb

Sign + or -	Eight	data	on	weight	includ- ing	the	deci- mal	point	Space	Unit	weight	<CR>	<LF>
-----------------------	--------------	-------------	-----------	---------------	------------------------	------------	----------------------	--------------	--------------	-------------	---------------	-------------------	-------------------

Examples:

+ 43.0913 g<CR><LF>

+1.519500 oz<CR><LF>

+215.3465 ct<CR><LF>

+ 0.094920 lb<CR><LF>

In part counting mode

Seven		data		on		parts	Space	p	c	s	<CR>	<LF>
--------------	--	-------------	--	-----------	--	--------------	--------------	----------	----------	----------	-------------------	-------------------

Example:

47 pcs <CR><LF>

In percentage mode:

Eight	data	on	percentages	includ- ing	the	decimal	point	Space	%	<CR>	<LF>
--------------	-------------	-----------	--------------------	------------------------	------------	----------------	--------------	--------------	----------	-------------------	-------------------

469.87 %<CR><LF>

12. MAINTENANCE

CONSERVATION IN GOOD WORKING ORDER, WASTE TREATMENT

Before undertaking any maintenance, cleaning or repair work, disconnect the device from the power supply.

12.1 CLEANING

Do not use aggressive agents (solvents, etc.). Clean with a cloth moistened with soapy liquid. Liquids must not get inside the device. Dry with a soft dry cloth.

Loose residues of the samples/dust can be removed with a brush or a hand-held vacuum cleaner.

If any material is spilled, remove it immediately.

12.2 MAINTENANCE, CONSERVATION IN GOOD WORKING ORDER

- The device may only be operated and maintained by trained and authorized Baxtran personnel.
- Before opening the device, it must be disconnected from the power supply.

12.3 WASTE TREATMENT

The packaging and the device must be recycled in accordance with the national or regional laws in force at the device's place of use.

13 HELP IN CASE OF MINOR BREAKDOWNS

In the event of alterations in the operation of the scale's program, turn it off and disconnect it from the power supply for a short time. Afterwards, the weighing process can be started again.

HELP	
Breakdown	Possible causes
The mass indicator does not work	• The scale is switched off. • The connection to the power supply is missing (unconnected/damaged power cable). • Power failure in the power grid.
The weight display is constantly fluctuating.	• Air currents/drafts • Vibration from the table/ground • The scale's plate is in contact with foreign bodies. • Due to electromagnetic fields/electrostatic charges (choose another location for installing the scale/if possible, turn off the device that is causing the disturbance)
The weighing result is obviously wrong	• The scale's display has not been reset to zero. • Incorrect calibration. • Problems with levelling the scale. • There are strong variations in temperature. • The preparation time was not respected. • Due to electromagnetic fields/electrostatic charges (choose another location for installing the scale/if possible, turn off the device that is causing the disturbance)

13.1 ERROR MESSAGES

message	cause	solution
E	The limit of the weighing range has been exceeded, excessive load.	Unload the scale.
-E	Weighing range not reached	e.g. the plate is missing
CAL-no	Settings' error	

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